
ORALB

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Oral-B®

powered by **BRAUN**

PRO

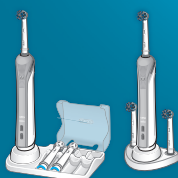Charger Type 3757
Handle Type 3766

2769271











- (UK)** 0 800 731 1792
- (IE)** 1 800 509 448
- (FR)** 0 800 944 803
(Service et appel gratuits)
- (BE)** 0 800 14 592
- (LU)** 800 21172
- (ES)** 900 814 208
- (PT)** 808 20 00 33
- (IT)** 800 124 600
- (NL)** 0800-4 45 53 88
- (PL)** 801 127 286
801 1 BRAUN
- (TR)** 0850 220 0911
- (GR)** 800 801 3457
- (ZA)** **Consumer Care Line 0860 112 188**
(Sharecall charged at local rates)

Internet:

www.oralb.com
www.braun.com
www.service.braun.com

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Charger type 3757
Handle type 3766

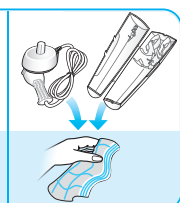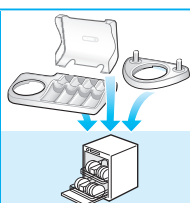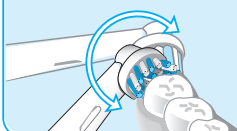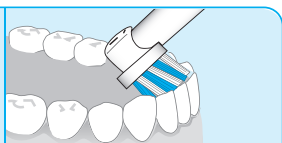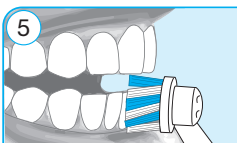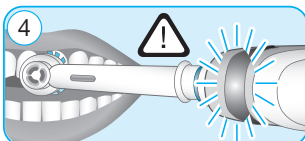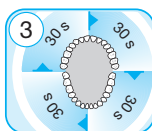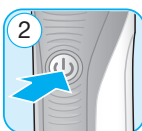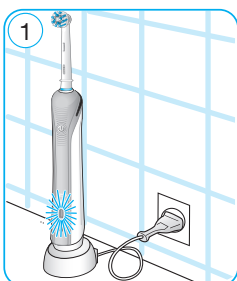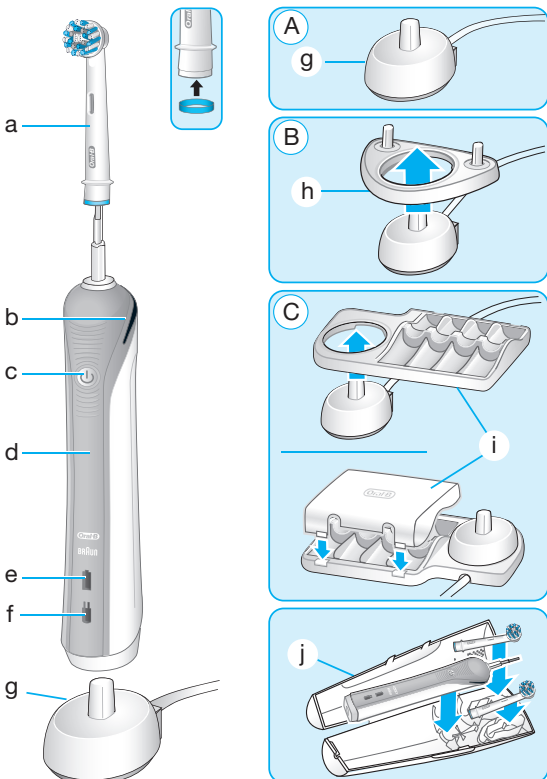English

Your **Oral-B®** toothbrush has been carefully designed to offer you and your family a unique brushing experience that is both safe and effective.

IMPORTANT

- Periodically check the entire product/cord for damage. A damaged or non-functioning unit should no longer be used. If the product/cord is damaged, take it to an Oral-B Service Centre. Do not modify or repair the product. This may cause fire, electric shock or injury.
- Usage by children under age 3 is not recommended. Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Cleaning and maintenance shall not be performed by children.
- Children shall not play with the appliance.
- Use this product only for its intended use

pressure. In addition the pulsation of the brush head will stop and the oscillation movement of the brush head will be reduced (during «Daily Clean», and «Sensitive» mode) (picture 4). Periodically check the operation of the pressure sensor by pressing moderately on the brush head during use.

Brush heads

Most Oral-B brush heads feature light blue INDICATOR® bristles to help you monitor brush head replacement need. With thorough brushing, twice a day for two minutes, the blue color will fade halfway approximately within 3 months, indicating the need to exchange your brush head. If the bristles splay before the color recedes, you may exert too much pressure on teeth and gums.

We do not recommend using the Oral-B «FlossAction» or Oral-B «3D White» brush head with braces. You may use the Oral-B «Ortho» brush head, specifically designed to clean around braces and wires.

Cleaning recommendations

After brushing, rinse your brush head under running water with the handle switched on. Switch off and remove the brush head. Clean handle and brush head separately; wipe them dry before reassembling. Unplug the charging unit before cleaning. Brush head compartment/ protective cover (i) and brush head holder (h) are dishwasher-safe. The charger (g) and travel case (j) should be cleaned with a damp cloth only (picture 6).

Make sure to store only cleaned and dried handle/brush heads in the travel case (j).

Subject to change without notice.

Environmental Notice

Product contains batteries and/or recyclable electric waste. For environment protection do not dispose in household waste, but for recycling take to electric waste collection points provided in your country.

Guarantee

We grant 2 years limited guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Service Centre.

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